



## Crutch clamp for 20mm diameter crutches



Kalle Fornication

[VIEW IN BROWSER](#)

updated 18. 2. 2024 | published 18. 2. 2024

### Summary

Click-in your crutches to prevent them from tipping over.


[Healthcare](#) > [Home Medical Tools](#)

Tags: [crutch](#) [crutchholder](#)

Click-in clamp for 20mm crutches on one side, screw fixation on the other side

You will need two threaded inserts M4 and two countersunk screws M4x20. No support needed. Print in PETG, it needs some elasticity. Pressing in a crutch on the screw side is quite hard, but you need to do that only once. Print the clamp, insert the threaded inserts, press in one crutch in the screw side, attach the distance piece and screw together. The other crutch can be clicked-in.

# Model files

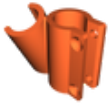

**clamp.stl**

☐ print in PETG



**distance\_piece.stl**

## License

This work is licensed under a  
[Creative Commons \(4.0 International License\)](https://creativecommons.org/licenses/by-nc/4.0/)



**Attribution-NonCommercial**

- 
- ✗ | Sharing without ATTRIBUTION
  - ✓ | Remix Culture allowed
  - ✗ | Commercial Use
  - ✗ | Free Cultural Works
  - ✗ | Meets Open Definition